

Etac Light cutlery

Award-winning ergonomically designed aids for daily living



Light is a lightweight cutlery series with functional design which facilitates the use and eating when strength and mobility in hands and arms is restricted. The cutlery is available in various models and is dishwasher-safe.

Cutlery with thick handles



The cutlery is lightweight, gentle on the joints and relieves strain on the fingers. The handles are oval to sit securely in the hand and allows various grips. The spoon has a slightly longer and thinner handle to make a grip in the hollow between the thumb and index finger comfortable. To make cutting easier the knife handle has a slope for the thumb to rest and a sharp blade.

Description Code	Dimen	sions	Weight	Or	der
Light Fork, Thick Handle	180 mn	า	27 g	80	402002
Light Spoon, Thick Handle	220 mn	ı	37 g	80	402004
Light Knife, Thick Handle	210 mn	ı	23 g	80	402006
Light Dessert Spoon, Thick Handle	180 mn	ı	30 g	80	402007
This cutlery relieves the strain on the joints in the hand and fingers.		Materia	I		Care
		Utensil parts: Stainless steel Handle: ABS plastic		eel	Dishwasher-safe

Spoon, angled

The lightweight angled spoons are designed for people with impaired arm and hand mobility. The handles are oval to sit in the hand and allow a secure grip when eating.






Cutlery with thin handles



This lightweight cutlery requires a minimum of effort to use. It has long thin handles to allow a comfortable grip in the hollow between the thumb and index finger. The long handles reduce the need to lift the arm as high as with normal length cutlery. The underside is gently rounded to prevent point pressure.

This cutlery is light and require a minimum of effort to use.

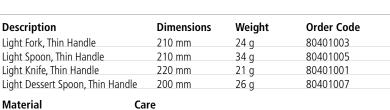

Utensil parts: Stainless steel Dishwasher-safe Handle: ABS plastic





Combination cutlery

This lightweight cutlery consists of a knife/fork and knife/spoon and is designed for use with one hand only. The handle permits a secure grip, with its rounded underside and flat upper side. The cutting edge is on the underside of the blade to prevent injury to the corners

of the mouth. The knife/fork is available in either righthanded and left-handed versions, while the knife/ spoon is suitable for both right-handed and lefthanded use.

Description	Dimensions	Weight	Order Code	
Light Knife/Spoon	200 mm	35 g	80403001	
Light Knife/Fork, Right Hand	180 mm	25 g	80403002	
Light Knife/Fork, Left Hand	180 mm	25 g	80403003	
5		5		

Material Utensil parts: Stainless steel Handle: ABS plastic **Care** Dishwasher-safe The combination cutlery provides a good grip for both cutting and picking up the food.





pursue t

"Everyone should be able to live a free and independent life and pursue their dreams regardless of any physical circumstances."