



# JOYCEeasy

*... is the intelligent choice.  
Nasal or Full Face. Simply easy.*

New all over: The delicate yet stable shape and fastening and setting elements on our new ventilation mask JOYCEeasy are designed for the simplest handling ever.

#### Intuitively simple use

In the nasal or full face variants, JOYCEeasy is easy to use and comfortable to wear. The intuitive use starts with adjustments made to the mask while it is on the patient.

#### Improved therapy compliance

JOYCEeasy also improves therapy compliance and the effectiveness of treatment. Several product tests have shown that we successfully optimized mask fit while minimizing the risk of leakage and pressure points.

#### Feelings of comfort

The WEINMANN SilkTec coating on JOYCEeasy makes the mask particularly comfortable. The special silky-soft surface on the mask cushion and forehead cushion provides greater wearing comfort than conventional silicone products.

#### Simple use

- Put on and take off mask with one hand
- Headgear clips smoothly click into place **4**
- Improved, even simpler fitting
- Hook-and-loop tabs and inner and outer sides of headgear are color-coded **5**

#### Simply comfortable



- Forehead cushion and mask cushion coated with WEINMANN SilkTec increase wearing comfort **1**
- The ball-and-socket joint lets patients move in every direction **6**
- Quiet exhalation system
- Adjustable cross-band keeps mask from slipping out of place **8**
- Large and soft around nasal root **1**

#### Attention to detail

- Enhanced fit reduces risk of leakage or pressure points to a minimum.
- Tube anchoring strap decreases leverage on tube
- Polyamide makes mask extremely robust
- Quick-release cord (only Full Face) **9**



**JOYCEeasy**  
Full face and nasal masks

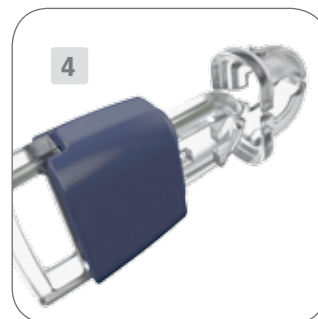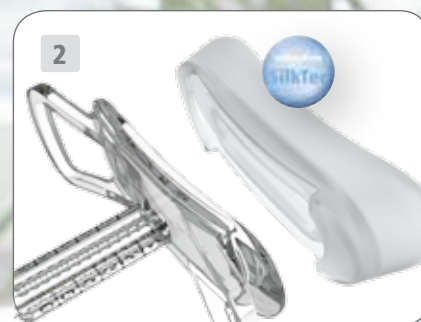

- 1 Mask cushion
- 2 Forehead cushion
- 3 Forehead support adjustment
- 4 Headgear clips
- 5 Hook-and-loop tab
- 6 Ball-and-socket joint
- 7 Exhalation system
- 8 Cross-band
- 9 Quick-release cord (for Full Face only)
- 10 Pressure measurement tube



**JOYCEeasy**  
Nasal mask



# JOYCEeasy Nasal








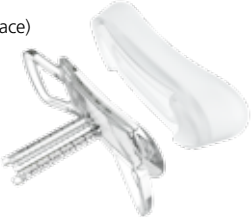
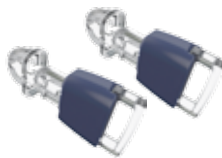





# Full Face

JOYCEeasy Full Face  
Full face mask




## Patient Interface / JOYCEeasy

WM Nr.			
Product name	WM Nr. Size S	WM Nr. Size M	WM Nr. Size L
JOYCEeasy Nasal	WM 25510	WM 25520	WM 25530
JOYCEeasy Full Face	WM 25910	WM 25920	WM 25930

Replacement parts			
<b>Headgear</b> incl. headgear clips WM 25505 (Nasal and Full Face)		<b>Forehead cushion,</b> silicone, WEINMANN SilkTec coating WM 26801 (Nasal and Full Face)	
<b>Tube anchoring strap</b> WM 26774 (Nasal and Full Face)		<b>Forehead support</b> incl. silicone forehead cushion, WEINMANN SilkTec coating WM 15797 (Nasal and Full Face)	
<b>Headgear clips</b> 2 pieces WM 26848 (Nasal and Full Face)		<b>Mask cushion</b> Silicone, WEINMANN SilkTec coating Size S WM 26803 (Nasal) WM 26853 (Full Face) Size M WM 26813 (Nasal) WM 26863 (Full Face) Size L WM 26823 (Nasal) WM 26873 (Full Face)	
<b>Quick-release cord</b> WM 26890 (Full Face)		<b>Ball-and-socket set, vented nasal</b> retaining ring, elbow, rotating sleeve WM 15798	
		<b>Ball-and-socket set, vented full face</b> retaining ring, elbow, anti-asphyxia valve, rotating sleeve WM 15897	
		<b>Carrying bag</b> WM 26295 (Nasal and Full Face)	

Our complete selection of therapy solutions, accessories and mask systems can be found at: [weinmann.de](http://weinmann.de)

Technical Data JOYCEeasy			CE 0197			
	Nasal mask	Full face mask	Nasal masks		Full face masks	
<b>Product class as per Directive 93/42/EWG</b>	II a					
<b>Dimensions (W x H x D)</b>	about 80 x 130 x 90 mm <sup>1</sup>	about 120 x 150 x 110 mm <sup>1</sup>	Inspiration at 50 l/min: Expiration at 50 l/min:		0.5 hPa 0.5 hPa	
<b>Weight</b>	about 72 g <sup>1</sup>	about 103 g <sup>1</sup>	<b>Switching pressure for Anti-asphyxia valve</b>		■ Open: ■ Close:	
<b>Therapy pressure range</b>	4 hPa to 30 hPa		-		≤ 1 hPa ≤ 2 hPa	
<b>Tube connection</b> Tapered connection as per EN ISO 5356-1	Ø 22 mm (male)		<b>Quoted two-figure noise emission value as per ISO 4871</b>		■ Sound pressure level: ■ Sound power level: ■ Uncertainty factor:	
<b>Pressure measurement tube</b>	Ø 4 mm				18 dB(A) 23 dB(A) 3 dB(A)	
<b>Temperature range</b>	■ Operation: +5° C to +40 °C ■ Storage: -20° C to +70 °C		<b>Service life</b>		up to 12 months <sup>2</sup>	
<b>Flow resistance</b>	■ at 50 l/min: ■ at 100 l/min:	0.16 hPa 0.64 hPa	<b>Standards applied</b>		EN ISO 17510-2:2009	
		0.3 hPa 0.53 hPa				

<sup>1</sup> Depending on size and design

<sup>2</sup> Materials used in the manufacture of masks age when exposed to aggressive cleaning products. Therefore, it may sometimes be necessary to replace mask parts earlier.

[weinmann.de](http://weinmann.de) | [weinmann.ch](http://weinmann.ch) | [weinmann-france.fr](http://weinmann-france.fr) | [weinmann.ru](http://weinmann.ru) | [weinmann.cn](http://weinmann.cn) | [weinmann.co.il](http://weinmann.co.il)

home care



# JOYCEeasy

Easy to Fit, Easy to Use