

made for you!



JOYCEone



JOYCEane

Your number 1 choice for a perfect fit in no time.

Time-consuming settings are history. Today we offer you JOYCEone, the ventilation mask in one size that fits nearly everyone. This mask, the result of years of experience and mask expertise, is available in a nasal or full face model.

Hand-in-hand: comfort and effective therapy

The forehead support's spring action automatically ensures that the proper distance is maintained for patient comfort and effective therapy. The single layer mask cushion is just as cleverly designed. It finds the ideal position on its own, regardless of the size and shape of the patient's nose and face.

Design is in the details

Your patients will love JOYCEone! This mask features more than a perfect fit and wearing comfort. The unmistakable color-coding of headgear and clips and the convenient recessed grip help your patient to assemble JOYCEone quickly and correctly right from the start.

Modern Efficiency

Short and sweet: the intelligent JOYCEone helps you to save time without compromising on patient care. JOYCEone – made for you!

The advantages of JOYCEone


- Saves time: automatic adjustment of forehead support 3 and mask cushion 1
- Unmistakable assembly: color-coding on headgear and clips 4
- Economical: one size fits nearly all patients
- Convenient: effortless cleaning
- Pleasant: quiet exhalation system
- Stable: very durable materials


JOYCEone Nasal

Mask fit on different nose sizes

JOYCEone Full Face Mask fit on range of nose

and face sizes

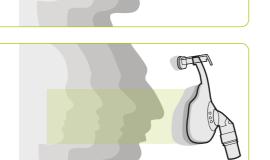

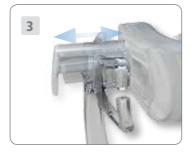

- 1 Mask cushion
- 2 Forehead cushion
- 3 Forehead support
- 4 Headgear clips
- 5 Hook-and-loop tab
- 6 Pressure measurement port
- 7 Quick-release cord (only JOYCEone Full Face)









Full Face



Patient Interface / JOYCEone



You'll find more information about our therapy solutions, accessories and mask systems at weinmann-medical.com

		JOYCEone Nasal	JOYCEone Full Face		JOYCEone Nasal	JOYCEone Full Face
	Product class according to 93/42/EEC lla	lla		Flow resistance ■ at 50 l/min ■ at 100 l/min	0.05 hPa 0.28 hPa	
	Dimensions (W x H x D)	ca. 96 x 115 x 83 mm	ca. 120 x 150 x 100 mm	Anti-asphyxia valve	0.2	
	Weight	87 g	100 g	flow resistance Inspiration at 50 l/min: Exhalation at 50 l/min:	-	0.5 hPa 0.5 hPa
	Therapy pressure range	4 bis 30 hPa		Switching pressure for anti-asphyxia valve © Open: © Close:		
	Hose connection: Cone according to EN ISO 5356-1	Ø 22 mm (male)				1 hPa 2 hPa
	Pressure measurement port	Ø 4 mm		Quoted two-figure noise emission value according to ISO 4871:	15 dB(A) 23 dB(A) 3 dB(A)	
,	Width of headgear clips	24 mm		Sound pressure level: Sound power level: Uncertainty factor:		
	Temperature range ■ Operation:	+5 °C to +40 °C -20 °C to +70 °C		Service life	up to 12 months ¹	
1	Storage:			Standards applied	EN ISO 17510-2:2009	

¹ Materials used in the manufacture of masks age when exposed to aggressive cleaning products. Therefore, it may sometimes be necessary to replace mask parts earlier.

Product name	WM no.
JOYCEone Nasal	WM 25260
JOYCEone Full Face	WM 25290

CECTIFIED QUALITY Management System 9342/EEC Arnex II (EN ISO 13485)

weinmann-medical.com I weinmann.de I weinmann.ch I weinmann.cn I weinmann.co.il I weinmann-france.fr I weinmann.se



JOYCEane

